
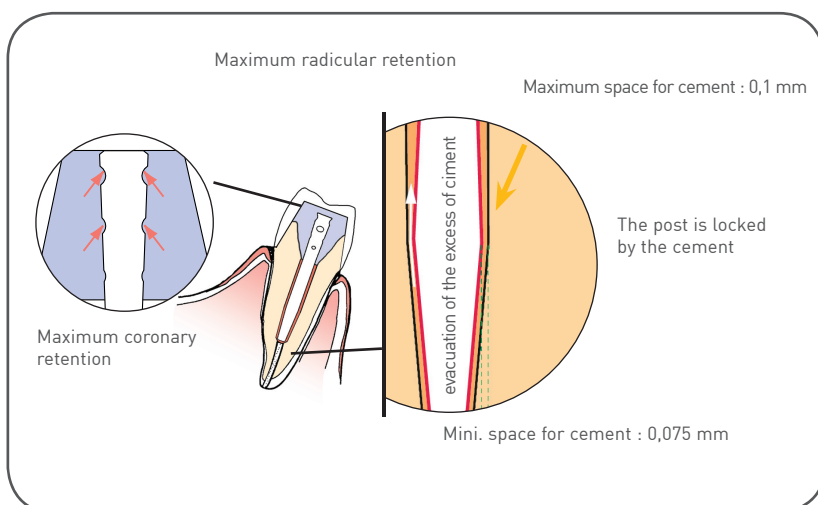


FIBIO® GLASSFIBER POSTS

Anthogyr has tested and selected one of the best glass fibers in terms of resistance, flexibility, adhesion properties and radio-opacity to produce its new Fibio® post. It can be used for both direct and indirect aesthetic restorations, without shadow under the ceramic.

- Even distribution of the stresses through the dentine, which significantly reduces the risk of fracture.
- Strong adhesion as the glassfiber is chemically close to the adhesion materials so as to get a global homogeneous restoration.
- Unique shape offering maximum retention in the coronary part as in the radicular part.

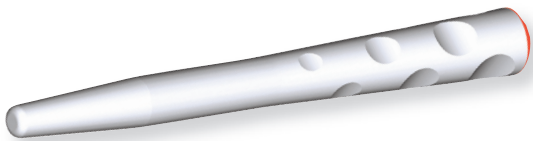

INTRODUCTION KIT

Réf. 9800



COMPOSITION

- Alkali / acid-resistant, unidirectional (longitudinal), zircon-reinforced silica fibres embedded in bisphenol A epoxy resin.
- Fiber : 43 %
- Zircon : 17 %
- Resin : 40 %


ELASTICITY

- Axial : 42,5 GPa
- Transversal : 15 GPa


CLEVER DISPENSER

- Selective opening and intelligent storage.
- Ergonomic dispenser : can be used with one hand.
- Easy selection of the size thanks to the calibrating film.


MODULE D'ELASTICITE EN TRACTION (MPa)
ELASTICITY MODULUS IN TRACTION [MPa]

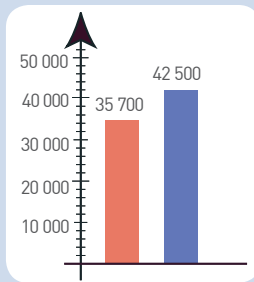

ESSAI DE FLEXION 3 POINTS (Flèche en mm)
BENDING STRESS TEST 3 POINTS (Arrow in mm)

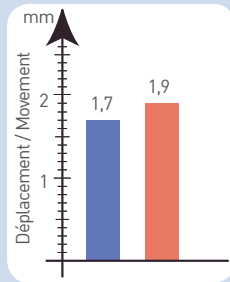

ESSAI DE COMPRESSION (Déformation en %)
COMPRESSION TEST (Deformation in %)

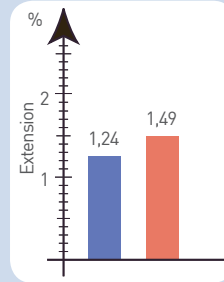

■ Fibres classiques
Classical fibres

■ Fibres renforcées
Reinforced fibres
(anthogyr)

	Foret ¹ Drills	Recharges Fibio ² Refills	Tenons d'empreinte ³ Impression posts
∅ 1,1	2891 T	FBC 1 lg 14	2881 IV
∅ 1,25	2892 T	FBC 2 lg 15	2882 IV
∅ 1,35	2893 T	FBC 3 lg 17	2883 IV
∅ 1,50	2895 T	FBC 4 lg 18	2885 IV



¹ Par 3 dans la taille souhaitée - Packs of 3 in the selected size.

² Par 10 ou par 100 de la même taille. Pour une recharge de 10, ajoutez «D» derrière la référence. Pour une recharge de 100, ajoutez «C» derrière la référence.
Packs of 10 or 100 in the selected size. For a refill of 10, add «D» at the end of the reference. For a refill of 100, add «C» at the end of the reference.

³ Par 20 de la même taille - Packs of 20 in the selected size.

	Ref. FBC 1	Ref. FBC 2	Ref. FBC 3	Ref. FBC 4	Ref. 2891	Ref. 2892	Ref. 2893	Ref. 2895	stop d'alésoir reamer stop	1 film de calibrage 1 calibrating film
Réf. 9800	x 5	x 5	x 5	x 5	x 1	x 1	x 1	x 1	x 20	•
Réf. 5608Q									x 50	
Réf. 2795										•